# *AccuProbe*

## **MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST**

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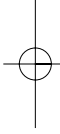
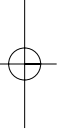
# AccuProbe

English

## MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST FOR EXPORT USE ONLY

(bioMérieux ref. 39004 / Gen-Probe Cat. No. 102840/2840)

### **INTENDED USE**

The ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST is a rapid DNA probe test which utilizes the technique of nucleic acid hybridization for the identification of *Mycobacterium intracellulare* isolated from culture.

### **SUMMARY AND EXPLANATION OF THE TEST**

*Mycobacterium intracellulare* (*M. intracellulare*) is a member of the *Mycobacterium avium* complex (*M. avium* complex) which consists of a number of organisms whose taxonomic relationships are both unclear and controversial, but whose pathogenicity in man is unquestioned (16). *M. intracellulare* has been shown to cause significant disease in immunocompromised patients (12). Treatment of this infection is difficult and the severity of the infection requires rapid diagnosis. Additionally in some laboratories the incidence of *M. avium* complex is equal to or greater than the incidence of *M. tuberculosis*.

Classical methods for identification of mycobacteria rely on staining specimens for acid fast bacilli followed by culture and subsequent biochemical testing. It can take as long as two months to speciate a *Mycobacterium* isolate using these standard culture methods (8).

*M. avium* complex is generally thought to consist of two species: *M. avium* and *M. intracellulare*. Phenotypically *M. avium* and *M. intracellulare* are virtually indistinguishable and biochemical tests are unable to differentiate between them.

High Performance Liquid Chromatography (HPLC) has been useful in the identification of *M. avium* and *M. intracellulare* (3), however, it is time consuming and not readily available to most clinical laboratories.

Serology has also been used for the differentiation of *M. intracellulare* and *M. avium* strains, using  $\alpha$ -antigen sera and may be useful in epidemiological studies. However, serotyping is not generally available and is of limited use in patient management. Currently there are 28 serovars generally accepted within the *M. avium* complex and at various times different serotypes have been assigned to the individual species *M. avium* and *M. intracellulare*. Historically serovars 1 through 3 were considered *M. avium*, while serovars 4 through 28 were considered *M. intracellulare* (15).

In several studies, Baess utilized DNA:DNA hybridization to clarify the taxonomic relationships between these two species. Based upon her analyses, she concluded that serovars 4, 5, 6, and 8, which at that time had been classified as *M. intracellulare*, actually belonged to the species *M. avium*. The status of serovar 9 was unclear (1, 2).

Saito, *et al*, also utilized DNA probe technology and has proposed reassigning the 28 serovars as follows:

serovars 1 through 6, 8 through 11, and 21 to *M. avium*; serovars 7, 12 through 20, and 25 to *M. intracellulare*, serovars 22 through 28 (except for serovar 25) were considered to be heterogeneous and could not be assigned to either species (13). In addition, some strains could not be serotyped, and a few agglutinated in more than one antisera.

Other reports utilizing DNA probes for speciating within the *Mycobacterium avium* complex have also been published, including the use of DNA probes for epidemiological studies and the geographical distribution of *M. avium* and *M. intracellulare* (5-7, 9, 11, 14).

The ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST identifies *M. intracellulare* isolated from culture in less than one hour of sample preparation. Identification is based upon the detection of specific ribosomal RNA sequences that are unique to *M. intracellulare*. The ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST offers a rapid, non-subjective and accurate means of identifying *M. intracellulare* isolated from culture.

### **PRINCIPLES OF THE PROCEDURE**

Nucleic acid hybridization tests are based on the ability of complementary nucleic acid strands to specifically align and associate to form stable double-stranded complexes (10). The ACCUPROBE SYSTEM uses a single-stranded DNA probe with a chemiluminescent label that is complementary to the ribosomal RNA of the target organism. After the ribosomal RNA is released from the organism, the labeled DNA probe combines with the target organism's ribosomal RNA to form a stable DNA:RNA hybrid. The Selection Reagent allows for the differentiation of non-hybridized and hybridized probe. The labeled DNA:RNA hybrids are measured in the GEN-PROBE luminometer. A positive result is a luminometer reading equal to or greater than the cut-off. A value below this cut-off is a negative result.

### **REAGENTS**

Reagents for the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST are provided in three separate reagent kits:

#### **ACCUPROBE MYCOBACTERIUM INTRACELLULARE PROBE KIT**

|   |                |
|---|----------------|
| <u>Probe Reagent</u> (P)<br><i>Mycobacterium intracellulare</i> | (4 x 5 tubes)  |
| <u>Lysing Tubes</u> (LT)<br>Glass beads and buffer              | (1 x 20 tubes) |

#### **ACCUPROBE CULTURE IDENTIFICATION REAGENT KIT**

|   |           |
|---|-----------|
| <u>Reagent 1</u> (Lysis Reagent) (1)<br>buffered solution containing 0.04% sodium azide | 1 x 10 mL |
| <u>Reagent 2</u> (Hybridization Buffer) (2)<br>buffered solution.                       | 1 x 10 mL |
| <u>Reagent 3</u> (Selection Reagent) (3)<br>buffered solution                           | 1 x 60 mL |

#### **GEN-PROBE DETECTION REAGENT KIT**

|  |            |
|--|------------|
| <u>Detection Reagent I</u> (RI)<br>0.1% hydrogen peroxide in 0.001 N nitric acid | 1 x 240 mL |
| <u>Detection Reagent II</u> (RII)<br>1 N sodium hydroxide                        | 1 x 240 mL |

## **WARNINGS AND PRECAUTIONS**

- A. For *in vitro* diagnostic use.
- B. Use universal precautions when performing this assay (4).
- C. Use only for the identification of *M. intracellulare* isolated from culture.
- D. Use only supplied or specified laboratory ware.
- E. Culture handling and all procedural steps through the heat inactivation step should be performed in a Class II Biological Safety Cabinet.
- F. Reagents in this kit contain sodium azide, which may react with lead or copper plumbing to form potentially explosive metal azides. Upon disposal of these reagents, always dilute the material with a large volume of water to prevent azide buildup in the plumbing.
- G. Avoid contact of Detection Reagents I and II with skin, eyes and mucous membranes. Wash with water if these reagents come into contact with skin. If spills of these reagents occur, dilute with water before wiping dry.

## **STORAGE AND HANDLING REQUIREMENTS**

Probe Reagent Tubes must be stored in the foil pouches at 2° to 8°C. The Probe Reagent Tubes are stable in the unopened pouches until the expiration date indicated. Once opened, the pouch should be resealed and the tubes should be used within two months and prior to the expiration date.

Other reagents used in the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST may be stored between 2° and 25°C and are stable until the expiration date indicated.

## **DO NOT FREEZE THE REAGENTS.**

## **SAMPLE COLLECTION AND PREPARATION**

The ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST is designed to determine the identity of *Mycobacterium intracellulare* isolated from culture.

- A. **Solid Media Method.** Growth from appropriate solid media, such as Lowenstein-Jensen slants or Middlebrook 7H10 or 7H11 plates, suggestive of *M. intracellulare* may be tested. Samples may be tested as soon as growth is visible and during the subsequent sixty days of incubation.
  - 1. Growth can be removed with a 1 mL disposable plastic loop, a wire loop, or a disposable plastic needle. Swabs should not be used due to the small volume of liquid in which the cells are subsequently resuspended.
  - 2. Avoid taking any of the solid media with the cells.
  - 3. The operator may elect to inoculate another culture plate at this time to confirm the purity of the isolate.
- B. **Broth Culture Method.** Growth in Middlebrook 7H9 broth with turbidity equivalent to or greater than a McFarland 1 Nephelometer Standard may be tested with the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST. Pipette a 100 µL sample from the well mixed broth suspension into the Lysing Tube as described below.

## **MATERIALS PROVIDED**

The ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST  
(bioMérieux ref. 39004 / Gen-Probe Cat. No. 102840/2840)

Probe Reagent (P)  
Lysing Tubes (LT)

### **20 Tests**

4 x 5 tubes  
1 x 20 tubes

## **MATERIALS REQUIRED BUT NOT PROVIDED**

1 µL plastic sterile inoculating loops, wire loops, or plastic needles for selecting colonies.  
Control culture strains  
Water bath or dry heat bath\* (60° ± 1°C)  
Water bath or dry heat bath\* (95° ± 5°C)  
Micropipettes (100 µL, 300 µL)  
Re-pipettor (100 µL, 300 µL)  
Vortex mixer

\*Heating blocks in the dry heat bath should have wells that are correctly sized for 12 x 75 mm tubes. The use of GEN-PROBE dry heat baths is recommended.

## **AVAILABLE FROM YOUR GEN-PROBE DISTRIBUTOR:**

GEN-PROBE LEADER 50i Luminometer  
(bioMérieux ref. 39400 / Gen-Probe Cat. No. 103100i/3100i)  
ACCUPROBE CULTURE IDENTIFICATION REAGENT KIT  
(bioMérieux ref. 39305 / Gen-Probe Cat. No. 102800/2800)  
GEN-PROBE DETECTION REAGENT KIT  
(bioMérieux ref. 39300 / Gen-Probe Cat. No. 201791/1791)  
Dry Heat Bath (60° ± 1°C)  
(bioMérieux ref. 39406 / Gen-Probe Cat. No. 3397)  
Dry Heat Bath (95° ± 1°C)  
(bioMérieux ref. 39407 / Gen-Probe Cat. No. 3398)  
Twin Dry Heat Bath (60°/95° ± 1°C)  
(bioMérieux ref. 39408 / Gen-Probe Cat. No. 3399)  
GEN-PROBE Sonicator  
(bioMérieux ref. 39409 / Gen-Probe Cat. No. 901104/T460)  
GEN-PROBE Sonicator Rack  
(bioMérieux ref. 39313 / Gen-Probe Cat. No. 104027/4027)

## **TEST PROCEDURE**

### **A. EQUIPMENT PREPARATION**

1. For optimal transfer of sonic energy, water must be thoroughly degassed according to the following procedure:
  - a. Add enough hot water to fill the sonicator bath to within 1 cm of the top of the tank.
  - b. Run the sonicator for 15 minutes to thoroughly degas the water.
2. Adjust one heating block or water bath to 60° ± 1°C and another heating block or water bath to 95° ± 5°C.
3. Prepare the GEN-PROBE luminometer for operation. Make sure there is sufficient volume of Detection Reagents I and II to complete the tests.

## B. CONTROLS

1. Positive and negative control strains should be tested routinely in each laboratory according to local regulations. A culture of *M. intracellulare* (e.g., American Type Culture Collection, ATCC #13950) may be used as the positive control while a culture of *M. avium* (e.g., ATCC #25291) may be used as the negative control.

## C. SAMPLE PREPARATION

1. Label a sufficient number of Lysing Reagent Tubes to test the culture isolates and/or controls. Remove and retain the caps.
2. Pipette 100  $\mu\text{L}$  of Reagent 1 (Lysis Reagent) and 100  $\mu\text{L}$  of Reagent 2 (Hybridization Buffer) into all Lysing Reagent Tubes. **If broth cultures are to be tested, do not add Reagent 1 to the Lysing Reagent Tubes.**
3. Transfer the sample from the solid media or 100  $\mu\text{L}$  of a well mixed broth culture into the labeled Lysing Reagent Tubes as described in the SAMPLE COLLECTION AND PREPARATION Section. Twirl the loop or needle in the Reagent 1 and Reagent 2 diluent mixture to remove the cells if testing growth from solid media.
4. Recap the Lysing Reagent Tubes and briefly VORTEX.

## D. SAMPLE LYSIS

1. Push the Lysing Reagent Tubes through the Sonicator Rack so that the reaction mixture in the bottom of the tube is submerged but the caps are above the water. Place Sonicator Rack on water bath sonicator. **DO NOT ALLOW THE TUBES TO TOUCH THE BOTTOM OR SIDES OF THE SONICATOR.**
2. Sonicate for 15 minutes.
3. Place the Lysing Reagent Tubes containing the sonicated organisms in a heating block or water bath for 10 minutes at  $95^{\circ} \pm 5^{\circ}\text{C}$ .
4. Carefully remove the Lysing Reagent Tubes from the heating block or water bath.

## E. HYBRIDIZATION

1. Open the foil pouch by cutting evenly across the top of the pouch. Remove enough Probe Reagent Tubes to test the culture isolates and/or controls. Reseal the pouch by folding the opened edge over several times and securing with adhesive tape or a clip. **Leave the desiccant pillow in the pouch.**
2. Label a sufficient number of Probe Reagent Tubes to test the culture isolates and/or controls. Remove and retain the caps.
3. Pipette 100  $\mu\text{L}$  of the lysed specimens from the Lysing Reagent Tubes into the corresponding Probe Reagent Tubes.
4. Recap the Probe Reagent Tubes and incubate for 15 minutes at  $60^{\circ} \pm 1^{\circ}\text{C}$  in a water bath or heating block.

## F. SELECTION

1. Remove the Probe Reagent Tubes from the water bath or heating block. Remove and retain the caps. Pipette 300  $\mu\text{L}$  of Reagent 3 (Selection Reagent) into each tube. Recap the tubes and VORTEX them to mix completely.

2. Incubate the Probe Reagent Tubes for 5 minutes at  $60^{\circ} \pm 1^{\circ}\text{C}$  in a water bath or heating block.
3. Remove the Probe Reagent Tubes from the water bath or heating block and leave them at room temperature for at least 5 minutes. Remove and discard the caps. **Read the results in the luminometer within 1 hour.**

#### G. DETECTION

1. Select the appropriate protocol from the menu of the luminometer software.
2. Using a damp tissue or paper towel, wipe each tube to ensure that no residue is present on the outside of the tube and insert the tube into the luminometer according to the instrument directions.
3. When the analysis is complete, remove the tube(s) from the luminometer.

#### **PROCEDURAL NOTES**

- A. **REAGENTS:** Reagent 2 (Hybridization Buffer) may precipitate. Warming and mixing the solution at  $35^{\circ}$  to  $60^{\circ}\text{C}$  will dissolve the precipitates.
- B. **TEMPERATURE:** The Hybridization and Selection reactions are temperature dependent. Therefore, it is imperative that the water bath or heating block is maintained within the specified temperature range.
- C. **TIME:** The Hybridization and Selection reactions are time dependent. Hybridize at least 15 minutes but no more than 20 minutes. Incubate the Probe Reagent Tubes during the SELECTION Step for at least 5 minutes but no more than 6 minutes.
- D. **WATER BATH:** The level of water in the water bath should be maintained to ensure that the Lysing Reagent Tubes are submerged up to, but not above, the level of the sealing ring. It should also be ensured that the entire liquid reaction volume in the Probe Reagent Tubes is submerged.
- E. **VORTEXING:** It is critical to have a homogenous mixture during the SAMPLE PREPARATION and SELECTION steps, specifically after the addition of cells to Reagents 1 and 2 and after addition of Reagent 3.
- F. **TROUBLE-SHOOTING:**
  1. Elevated negative control values (*M. avium*, ATCC #25291) greater than 10,000 RLU (Relative Light Units) in the LEADER or 300 PLU (Photometric Light Units) in the ACCULDR (formerly PAL) can be caused by insufficient cell numbers, improper sonication, or by testing mixed or aged cultures. Because mixed cultures can occur, a portion of the growth may be streaked onto the appropriate agar medium and incubated to check for multiple colony types.
  2. Low positive control values (*M. intracellulare*, ATCC #13950) less than 30,000 RLU in the LEADER or 900 PLU in the ACCULDR (formerly PAL) can be caused by insufficient cell numbers, improper sonication, or by testing mixed or aged cultures. Because mixed cultures can occur, a portion of the growth may be streaked onto the appropriate agar medium and incubated to check for multiple colony types.

#### **RESULTS**

##### A. INTERPRETATION OF RESULTS

The results of the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST are based on the following cut-off values. Samples producing signals greater than or equal to these cut-off values are considered positive. Signals less than these cut-off values are considered negative. Results in repeat ranges should be repeated.

|               | <b>ACCULDR</b><br>(formerly PAL) | <b>LEADER</b>     |
|---------------|----------------------------------|-------------------|
| Cut-off Value | 900 PLU                          | 30,000 RLU        |
| Repeat range  | 600-899 PLU                      | 20,000-29,999 RLU |

#### B. QUALITY CONTROL AND ACCEPTABILITY OF RESULTS

Negative control (e.g., *M. avium*, ATCC #25291) and positive control (e.g., *M. intracellulare*, ATCC #13950) should satisfy the following values:

|                  | <b>ACCULDR</b><br>(formerly PAL) | <b>LEADER</b> |
|------------------|----------------------------------|---------------|
| Negative control | < 300 PLU                        | < 10,000 RLU  |
| Positive control | > 900 PLU                        | > 30,000 RLU  |

#### **LIMITATIONS**

This method has been tested using fresh growth from solid media and from broth cultures listed in the SAMPLE COLLECTION AND PREPARATION Section. The efficacy of this test has not been demonstrated on direct clinical specimens (e.g., urine, stool, or respiratory specimens).

Results from the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST should be interpreted in conjunction with other laboratory and clinical data available to the clinician.

#### **EXPECTED VALUES**

The ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST was compared to standard culture biochemical identification methods from two sites, site 1 and 2, using a total of 259 isolates: 109 isolates of the *M. intracellulare* and 147 isolates of 19 other *Mycobacterium* species. Standard culture identification is dependent on growth rate, colony morphology, microscopic examination, and a series of biochemical reactions. In addition, the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST was compared to High Performance Liquid Chromatography (HPLC) at site 3 using a total of 97 *Mycobacterium* strains. HPLC identified 31 isolates as *M. intracellulare*, 36 isolates as *M. avium*, and 36 as isolates representing 12 other *Mycobacterium* species. Using the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST all isolates were categorized as either positive (> 30,000 RLU) or negative (< 30,000 RLU). The range of observations for negative cultures was 266 to 3,405 RLU and 33,851 to 559,708 RLU for positive cultures. A comparison of these results to standard culture identification methods (Sites 1 and 2) and HPLC methods (Site 3) are shown below.

#### **ACCUPROBE / CULTURE AND HPLC IDENTIFICATION**

| <b>ACCUPROBE Culture</b> | <b>Pos Pos</b> | <b>Pos Neg</b> | <b>Neg Pos</b> | <b>Neg Neg</b> | <b>Sensitivity/<br/>Specificity</b> | <b>Percent Agreement</b> |
|--------------------------|----------------|----------------|----------------|----------------|-------------------------------------|--------------------------|
| Site 1                   | 28             | 0              | 0              | 50             | 100% / 100%                         | 100%                     |
| Site 2                   | 81             | 0              | 0              | 100            | 100% / 100%                         | 100%                     |
| Site 3                   | 31             | 0              | 0              | 66             | 100% / 100%                         | 100%                     |
| <b>Total</b>             | <b>140</b>     | <b>0</b>       | <b>0</b>       | <b>216</b>     | <b>100% / 100%</b>                  | <b>100%</b>              |

## **PERFORMANCE CHARACTERISTICS**

### **A. WITHIN-RUN PRECISION**

The within-run precision of the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST was calculated by assaying two concentrations of ribosomal RNA isolated from *M. intracellulare* using 10 replicates in a single assay.

#### *Mycobacterium intracellulare*

| <b>Sample</b>            | <b>A</b> | <b>B</b> |
|--------------------------|----------|----------|
| Number of Replicates     | 10       | 10       |
| Mean Response            | 39,179   | 71,587   |
| Standard Deviation       | 764      | 2,123    |
| Coefficient of Variation | 2.0%     | 3.0%     |

### **B. BETWEEN-RUN PRECISION**

The between-run precision was calculated by assaying the same two concentrations of *M. intracellulare* ribosomal RNA using single determinations in 12 consecutive runs.

#### *Mycobacterium intracellulare*

| <b>Sample</b>            | <b>A</b> | <b>B</b> |
|--------------------------|----------|----------|
| Number of Replicates     | 12       | 12       |
| Mean Response            | 37,541   | 72,189   |
| Standard Deviation       | 3,318    | 4,352    |
| Coefficient of Variation | 8.8%     | 6.0%     |

### **C. SPECIFICITY**

A total of 106 ATCC culture isolates were evaluated using the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST. These isolates represented a total of 91 species from 39 genera. Five isolates of *M. intracellulare*, 37 isolates of 25 other *Mycobacterium* species, and 64 isolates of 38 other genera representing a phylogenetic cross-section of organisms were evaluated using the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST. Only *M. intracellulare* isolates tested produced a positive result using the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST. Other *Mycobacterium* species and the representative phylogenetic cross-section species did not react using this kit.

### **D. RECOVERY**

*M. intracellulare* ribosomal RNA at concentrations ranging from  $5 \times 10^{-4}$   $\mu\text{g}$  and  $1 \times 10^{-1}$   $\mu\text{g}$  per test was assayed in the presence of 15 million cells of either *M. tuberculosis*, *M. simiae*, or *Nocardia asteroides*. No interference of *M. intracellulare* signal was observed and the other organisms present did not react using the ACCUPROBE MYCOBACTERIUM INTRACELLULARE CULTURE IDENTIFICATION TEST.

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

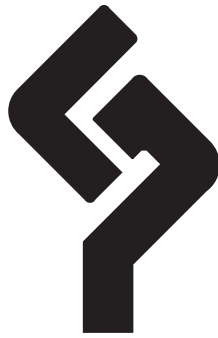
BIBLIOGRAPHY / LITERATUR / BIBLIOGRAPHIE / BIBLIOGRAFIA


1. **Baess, I.** 1983. Deoxyribonucleic acid relationships between different serovars of *M. avium*, *M. intracellulare* and *M. scrofulaceum*. ACTA Path. Microbiol. Scand. Sect. B. 91:201-203.
2. **Baess I.** 1979. Deoxyribonucleic acid relatedness among species of slowly-growing mycobacteria. ACTA Path. Microbiol. Scand. Sect. B. 87:221-226.
3. **Butler, W. R. and J. O. Kilburn.** 1988. Identification of major slowly growing pathogenic mycobacteria and *M. gordonae* by high performance liquid chromatography of their mycolic acids. J. Clin. Microbiol. 26:50-53.
4. **Centers for Disease Control.** 1988. United States Morbid. and Mortal. Weekly Rep. 37:377-382, 387-388.
5. **Drake, T. A., J.A. Hindler, G. W. Berlin and D. A. Bruckner.** 1987. Rapid identification of *M. avium* complex in culture using DNA probes. J. Clin. Microbiol. 25:1442-1445.
6. **Ellner, P. D., T. E. Kiehn, R. Cammarata and M. Hosmer.** Rapid detection and identification of pathogenic mycobacteria by combining radiometric and nucleic acid probe methods. J. Clin. Microbiol. 26:1349-1352.
7. **Gonzalez, R. and B. A. Hanna.** 1987. Evaluation of Gen-Probe DNA hybridization systems for the identification of *M. tuberculosis* and *M. avium-intracellulare*. Diagn. Microbiol. Infect. Dis. 8:69-77.
8. **Kent, P. T. and G. P. Kubica.** 1985. Public Health Mycobacteriology: A guide for the level III laboratory, U. S. Department of Public Health and Human Services. Public Health Service. Centers for Disease Control, Atlanta, GA.
9. **Kiehn, T. E. and F. F. Edwards.** 1987. Rapid identification using a specific DNA probe of *M. avium* complex from patients with acquired immunodeficiency syndrome. J. Clin. Microbiol. 25:1551-1552.
10. **Kohne, D.E., A.G. Steigerwalt and D.J. Brenner.** 1984. Nucleic acid probe specific for members of the genus Legionella, p. 107-108. In C. Thornsberry, et al. (ed.) Legionella: proceedings of the 2nd international symposium. American Society for Microbiology, Washington, D.C.
11. **Musial, C. E., L.S. Tice, L. Stockman and G. D. Roberts.** 1988. Identification of Mycobacteria from culture by using the Gen-Probe rapid diagnostic system for *M. avium* complex and *M. tuberculosis* complex. J. Clin. Microbiol. 26:2120-2123.
12. **Pitchenik, A. E., D. Fertel and A. B. Block.** 1988. Mycobacterial disease: epidemiology, diagnosis, treatment and prevention. Clin. Chest. Med. 9:425-441.
13. **Saito, H., H. Tomioka, K. Sato, H. Tasaka and D. Dawson.** 1990. Identification of various serovar strains of *M. avium* and *M. intracellulare*. J. Clin. Microbiol. 28:1694-1997.
14. **Saito, H., H. Tomioka, K. Sato, H. Tasaka, M. Tsukamura, F. Kuze and K. Asano.** 1989. Identification and partial characterization of *M. avium* and *M. intracellulare* by using DNA probes. J. Clin. Microbiol. 27:994-997.
15. **Schaefer, W. B.** 1965. Serologic identification and classification of the atypical mycobacteria by their agglutination. Am. Rev. Respir. Dis. 92(Supl.):85-93.
16. **Sommers, H. M. and R. C. Good.** 1985. M., p. 216-248. In E. H. Lennette et al. (ed.) Manual of Clinical Microbiology, 4th ed. American Society for Microbiology, Washington, D.C.

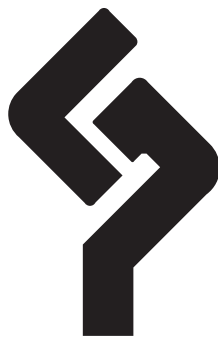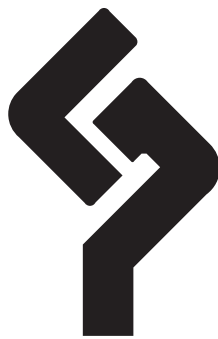**Gen-Probe Incorporated**  
**San Diego, CA 92121 (USA)**

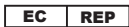Authorized Representative  
**EMERGO EUROPE**  
Molenstraat 15  
2513 BH The Hague  
The Netherlands

102901F-01 Rev. E 2011-03  
©1990 - 2011 Gen-Probe Incorporated