




DIACLONE Immunology Products


Anti-Human CD34 Class III PE Conjugated

Clone: B-K34

CATALOGUE N° :	854.702.010
CLONE :	B-K34
ISOTYPE :	Mouse IgG1
IMMUNISATION :	Human CD34 leukaemia cells
HYBRIDOMA :	Myeloma X63/AG.8653 x Balb/c spleen cells
SPECIFICITY :	Recognises the Gp105-120 antigen, a 105-120 kDa protein
WORKSHOP :	IV
QUANTITY :	100 tests
FORMAT :	Phosphate-buffered saline with 1% BSA and 0.1% sodium azide
WORKING DILUTION :	Reconstitute with 1 ml deionised water Use 10 µl to label 10 ⁶ cells or 100 µl of whole blood
STORAGE :	Stable at +2-8°C for six months after reconstitution DO NOT FREEZE
APPLICATION :	Suitable for flow cytometry

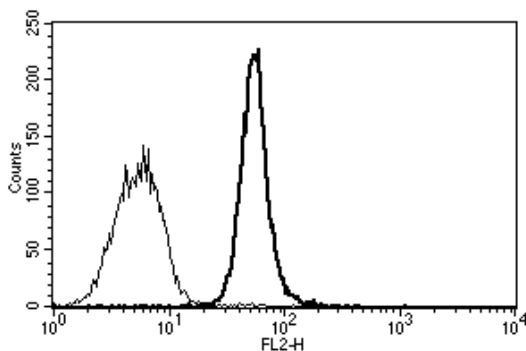

A typical staining pattern with the B-K34 monoclonal antibody of KG-1a cell line

 Sodium azide is a poisonous and hazardous substance, handle and dispose of properly

Version 1 - 05.01

For research use only

20 YEARS OF EXPERTISE

With over 20 years experience and extensive expertise we are committed to providing excellence in Monoclonal Antibody and Immunoassay development.

The expanding range of Diacclone Immunology products is specifically designed to advance research applications.

Our experience and expertise coupled to the diversity of our product range makes Diacclone a clear choice to

Fast Track Your Research...

DIACLONE PRODUCTS



Monoclonal Antibodies



ELISA



ELISpot



DIAPlex

GLOBAL DISTRIBUTION

For the USA please contact our distributor
Cell Sciences, Tel: 888-769-1246
email: info@cellsciences.com
For all other countries please visit our website
for a full list of our distributors.

